



Uncontrolled Copy when printed
CV-000918-TC, Revision 7



COVVI HAND & ASSOCIATED PRODUCTS

Technical Manual



**Anything Is
Possible**

Foreword

COVVI was founded in 2017 by a passionate and experienced team who wanted to continue to create leading multi-articulating device designs.

The human hand is amazing and the level of functionality which can be restored through prosthesis has historically been very low; but as engineers and designers we see this as a challenge to continuously improve our designs to meet the needs of users.

COVVI started with a blank sheet of paper and a mission to address some of the most prevalent deficiencies of current prosthetic hands. Fast forward to present day and we now have the best looking and most advanced multi-articulating hand on the market.

The COVVI Team



Contents

1.0	COVVI Product Range	6
2.0	Safety Precautions	8
2.1	Wrist Lamination	9
3.0	Environmental Conditions	10
4.0	The COVVI Hand	11
4.1	The Hand Box	11
4.2	Dimensions	13
4.3	COVVI Hand Overview	14
4.4	COVVI Wrist Variants	16
4.5	Technical Information	18
5.0	Fitting the COVVI Hand	19
5.1	Connecting/Disconnecting the COVVI Hand	19
5.2	E-Paper Screen	21
5.3	System Notifications	22
5.4	Glove Removal	23
6.0	Grip Patterns	24
6.1	Opposed Grips	24
6.2	Non-Opposed Grips	25
7.0	The COVVI Power Supply	26
7.1	What's Inside	27
7.2	Safety Precautions	28
7.3	COVPS-1600 Batteries	29
7.4	COVPS-2600 Batteries	29
7.5	Activating the System	30
7.6	Charge Indication	31
7.7	Fitting the COVVI Power Supply	32
7.8	System Specification	33
7.9	Battery Replacement	35

Contents

8.0	Electromagnetic Compatibility (EMC)	36
9.0	Component Compatibility	37
10.0	COVVI Go	40
10.1	COVVI Go	40
10.2	Registration	40
10.3	COVVI Portal	40
10.4	Hand Selection	41
10.5	Input Settings	41
10.6	Grip Tables	42
10.7	User Grips	43
10.8	Grip Triggers	43
10.9	Hand Control	45
10.10	Config Management	47
10.11	COVVI Touch	47
10.12	Games	48
10.13	Hand Statistics	48
10.14	System Information	48
10.15	Software Update	49
10.16	Remote Assist	49
11.0	Warranty	50
11.1	24-Month Manufacturer's Warranty	50
11.2	Optional Extended Warranty	50
11.3	Additional Information	50
12.0	Warranty Service and Repairs	51
12.1	Limited Liability	51
13.0	Device Problem Reporting	52
13.1	Returns	52
14.0	Disposal	53
15.0	Symbols Used	54
16.0	Declaration of Conformity	58


1.0 COVVI Product Range

Components for EQD (Flex and Straight Wrist Variants)

The Hand Box

Contains one unit of the COVVI Hand, housed in protective material and shipped in a custom-sized box. The device is available in either left or right orientations, and in three sizes: small, medium and large. The Carbon Black, White and Titan Grey colour options are available for all sizes, while the Rose Gold is an option for small and medium hands only.

The COVVI Hand shown here is the Flex Wrist version, it is also available with either a Straight Wrist or Friction Wrist variant.

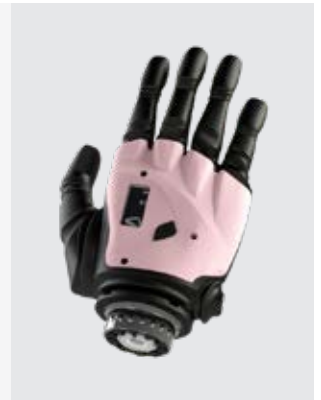

The Power Supply / Wall Charger Boxes

Spare Batteries and Wall Chargers can be ordered on their own, which will be shipped in their own custom boxes.

A full Battery Kit can also be ordered, which includes one power supply, a USB-C Wall Charger (with plug adaptors for different regions) and a USB-C Control Panel.

Please Note: Any orders for a Spare Battery or Battery Kit only contain either one of the 1600mAh or 2600mAh power supplies shown to the right. The Control Panel is also only included as a part of the Battery Kit.



COVPS-1600



COVPS-2600



Control Panel



Wall Charger
(with adaptors)

The Wrist Assembly Kit Package

(Also known as the 'COVVI Wrist Lamination Kit')

When ordering a Coaxial Core, Coupling Piece or one size of a Lamination Ring (40mm, 45mm, 50mm), these components will be shipped in a bio-degradable plastic pouch.

A Wrist Assembly Kit is available that bundles all three of these items together, where they will be packaged and shipped together.

Please Note: Multiple items from different product ranges can be ordered together, and they will all be securely packaged into one of our custom shipping boxes.



Coupling Piece



Coaxial Core



Lamination Ring

1.0 COVVI Product Range

The Glove Kit

Replacement Glove Kits can be ordered that contain a new COVVI Glove, Dorsal Cover and Retention Ring. Covers available in Carbon Black, White, Titan Grey and Rose Gold (Rose Gold is only available in small and medium). Please check the ordered kit is suitable for the size, orientation and type of COVVI Hand used.

Glove Kit: Compatible with the COVVI Hand.

Glove Kit 2: Comes with two Retention Ring halves that are compatible with the COVVI Hand - Straight Wrist and COVVI Hand - Friction Wrist variants.

Please Note: Fixings cannot be used across Glove Kit versions.



Components for Friction Wrist

COVVI Friction Wrist Switch Block

The Switch Block must be ordered alongside a COVVI Hand - Friction Wrist in order to connect the Hand to the auxiliary control and power components within the prosthesis socket.



COVVI Friction Wrist Lamination Ring

50mm Lamination Ring for the COVVI Hand - Friction Wrist.



COVVI Friction Wrist Installation Tool

An assistive tool that helps with the installation and removal of the COVVI Hand - Friction Wrist.



2.0 Safety Precautions



Please read the following safety precautions prior to fitting the COVVI Hand. Please note that the COVVI Hand is recommended with COVVI products.

In this Technical Manual, the term 'COVVI Hand' refers to all wrist variants. If any instructions, warnings, or information apply only to a specific wrist type, this will be clearly stated as applying to the COVVI Hand with Flex Wrist, Straight Wrist, or Friction Wrist.

The COVVI Hand should only be fitted by a Certified Prosthetist or Technician.

Certified Prosthetists and Technicians can remove the silicone glove from the hand and install a new silicone glove **ONLY** when training has been completed, and a certificate has been issued by a COVVI Representative.

NO other disassembly can be performed unless you have received specific training by COVVI and a certificate has been issued by a COVVI Representative.

COVVI Ltd. and COVVI USA INC have the right to void the warranty of all products that have any type of modification or damage caused by any unauthorised or untrained personnel.

- Do **NOT** use the COVVI Hand without the integrated silicone glove as it protects the hand from water ingress and other foreign bodies, which could damage the internal components and cause the hand to fail.
- The COVVI Hand has been tested and certified with an IP44 rating, classifying the degree of protection against solids and water of the hand to solid objects over 1mm in diameter and sprays and splashes from all directions. However, the wrist interface is not IP44 water or dust resistant.
- Do **NOT** submerge the COVVI Hand in water. Any damage caused by intentional submersion in water will not be covered under the warranty. If the glove is punctured, it will affect the IP44 rating of the hand.
- Do **NOT** attempt to use the COVVI Hand while the batteries are charging. When the batteries are charging, the power will automatically turn off. If for any reason the power does not turn off while charging, using the hand can be potentially unsafe.
- Do **NOT** use the COVVI Hand if a failure occurs, as it could cause damage to your patient or others.
- Do **NOT** use an additional glove with the COVVI Hand. Cosmetic gloves from other manufacturers restrict the movement of the hand, resulting in lower battery life, reduced grip strength and an increased internal temperature.
- The Essential Performance is to move the digits only when an input signal is applied. Should the fingers and/or thumb move without an input, the patient should stop using the COVVI Hand immediately and contact their clinician to get in touch with COVVI Customer Services.



This symbol is used in instances to warn about hazardous or dangerous materials and/or objects and will be used throughout this document.



2.0 Safety Precautions

- We do not recommend the use of carbon fibre materials in a fabricated socket unless it is grounded to the carbon fibre lamination.
- This hand uses semiconductors that can be damaged by electrostatic discharge (ESD).
- Individuals who are exposed to hazardous environments that contain flammable liquid or gas should **NOT** operate the COVVI Hand when in those environments.
- Do **NOT** subject the COVVI Hand to intentional excessive impacts. Any damage caused by intentional harm or neglect will not be covered under the warranty.
- Do **NOT** expose the COVVI Hand to an open flame.
- Do **NOT** use the COVVI Hand to operate a firearm.
- Do **NOT** expose the COVVI Hand to chemicals such as solvents, acids, alkalis, corrosive substances, detergents, and similar chemicals, as this could damage the components in the hand. If the COVVI Hand needs cleaning, mild soap and a soft damp cloth should be used. Abrasive cleaners and cleaning cloths will scratch the covers.
- Before operating a vehicle with the COVVI Hand, local regulations should be checked, and the patient's insurance company must be notified. We recommend that the prosthesis is fully charged. In the event of an accident occurring while operating a vehicle with a COVVI Hand or any other COVVI product, COVVI Ltd. and COVVI USA INC shall not be liable under any circumstances.
- The COVVI Hand is not suitable for extreme sports. These include, but are not limited to: base jumping, sky diving, downhill/endurance mountain biking, motocross, rock climbing, abseiling, bouldering and power lifting. Any sport that would involve the COVVI Hand getting wet must be avoided. These include, but are not limited to: diving, swimming, snorkelling, surfing, kayaking and canoeing.
- The COVVI Hand is **NOT** suitable to support body weight on crutches; we recommend that the hand is **NOT** used with crutches.
- Avoid remaining in the vicinity of visible strong magnetic or electrical interference, such as theft prevention systems.
- Maintain extra caution when handling sharp objects.
- Do **NOT** subject the hand to impact from sharp-edged or pointed objects.

COVVI Ltd. and COVVI USA INC have the right to void the warranty of all products that have any type of modification or damage caused by any unauthorised or untrained personnel.

2.1 Wrist Lamination

Ensure that the wrist is adequately laminated to the socket prior to use. Do not use the product if there are cracks in the lamination or if there is play between the socket and the Wrist Lamination Ring.

3.0 Environmental Conditions

The COVVI Hand is to be used exclusively for upper extremity device fittings. Based on the environmental conditions below, the patient should not use, transport, or store the hand outside of the boundaries listed below.


Maximum Temperature	104°F / 40°C
Minimum Temperature	-4°F / -20°C
Humidity	Maximum 80%, non-condensing
Pressure	101.3kPa / 14.7psi

4.0 The COVVI Hand

4.1 The Hand Box

COVVI Hand - Flex and Straight Wrist Packaging


To remove the COVVI Hand from its box, first slide foam piece 1 vertically out of the box, then slide the rest of the foam pieces along into the gap left by foam piece 1. Lift the hand from the box and remove foam piece 2 from the EQD and foam piece 5 from the fingers. To place the hand back in the box, follow the steps in reverse. Please take care with the orientation of the foam pieces.

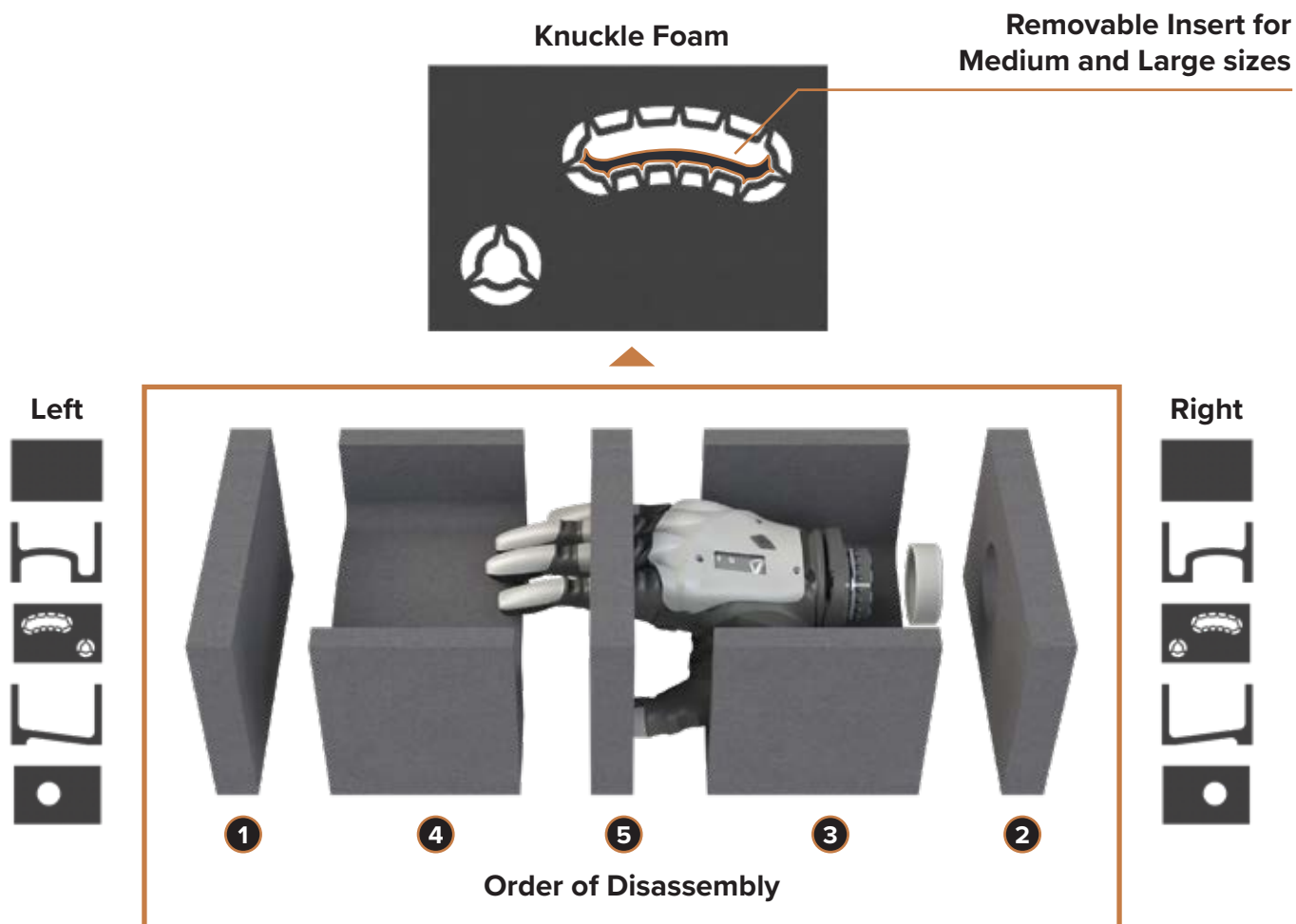

4.0 The COVVI Hand

COVVI Hand - Friction Wrist Packaging

To remove the COVVI Hand - Friction Wrist from its box, first slide foam piece 1 vertically out of the box, then slide the rest of the foam pieces along into the gap left by foam piece 1. Lift the hand from the box and remove foam piece 2 from the Friction Wrist Lamination, and foam piece 5 from the fingers. To place the hand back in the box, follow the steps in reverse. Please take care with the orientation of the foam pieces.

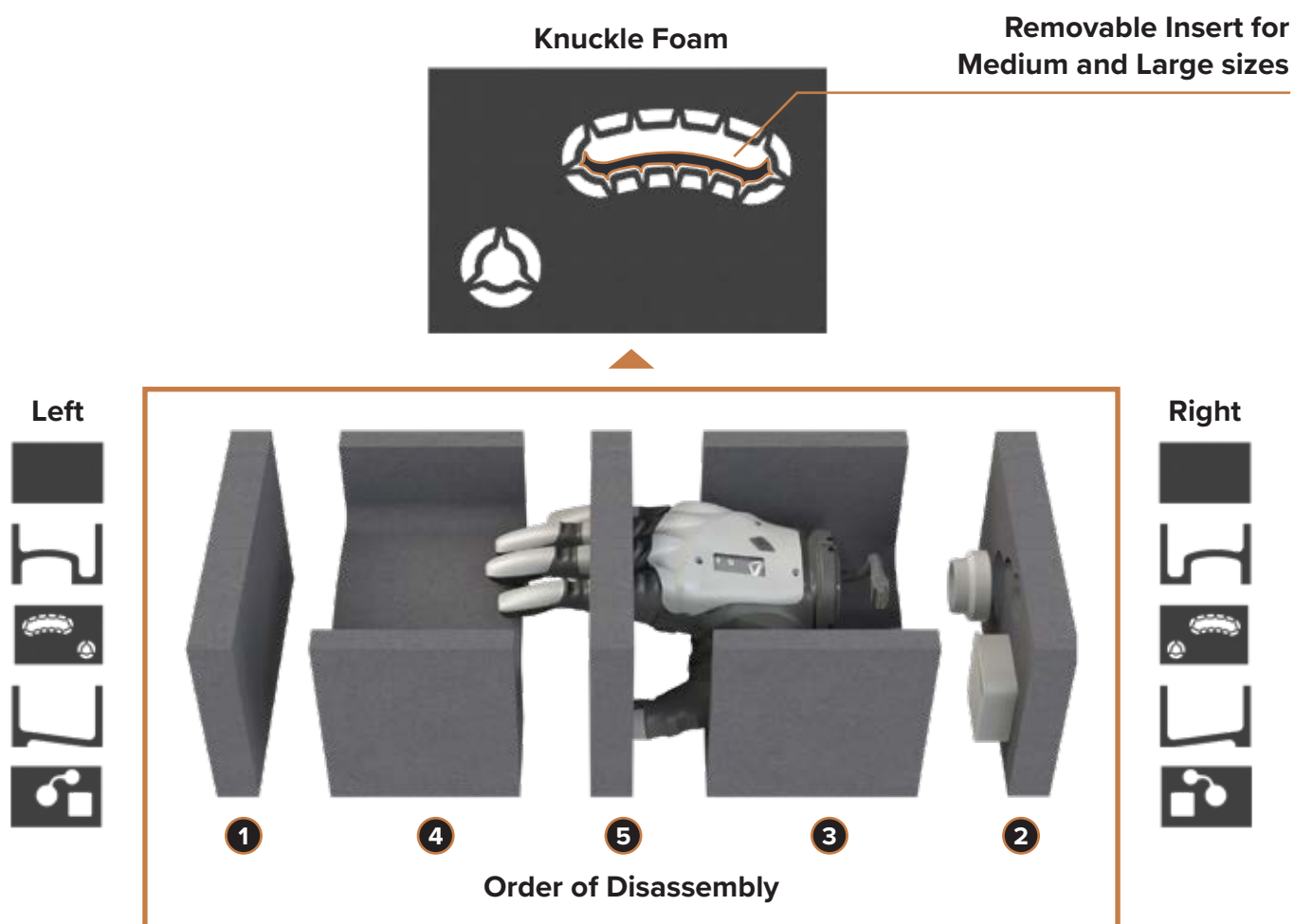

4.0 The COVVI Hand


4.2 Dimensions



Dimensions	Key	Small		Medium		Large	
		Inches	mm	Inches	mm	Inches	mm
Thumb Tip to Thumb Base	A	4.3	110	4.5	114	4.9	124
Middle Finger to Base (Flex Wrist)	B	6.6	168	7.1	179	7.4	189
Middle Finger to Base (Straight/Friction Wrist)	B	6.4	163	6.9	174	7.2	184
Max. Chassis Width	C	2.8	70	3	75	3.1	80
Palm Circumference	D	7.5	190	7.7	195	7.9	200
Maximum Opening Width - Tripod Grip	E	3.4	86	3.6	92	3.9	100

4.0 The COVVI Hand

4.3 COVVI Hand Overview

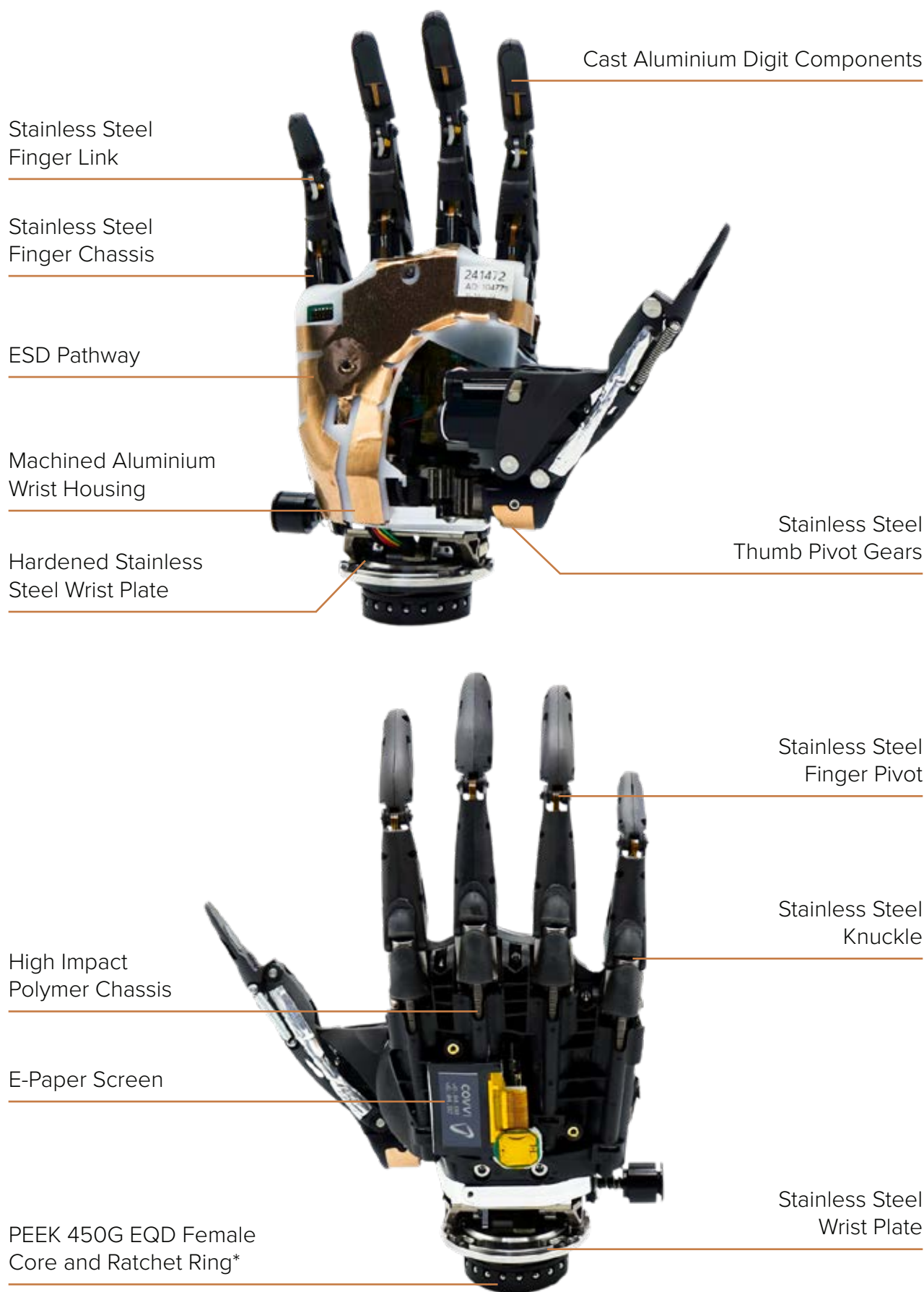

4.0 The COVVI Hand

Optimised Fingertip Geometry

Hardened 60A Shore Fingertips

E-Paper Screen

Dorsal Button

Wrist Lock Button**



Touchscreen Compatibility with the Index Finger

FSRs in all the Fingertips

Flexible Glove (IP44)

Thicker Palmar Region

Flex Wrist***

EQD - Electric Quick Disconnect*



Friction Wrist with Lamination Ring

*The EQD is not present if a Friction Wrist is installed

**The Wrist Locking Button is not operational if a Straight or Friction Wrist is installed

***The Flex Wrist is not present if a Straight or Friction Wrist is installed

4.0 The COVVI Hand

4.4 COVVI Wrist Variants

The COVVI Wrist Suite is comprised of three distinct variants of the COVVI Hand, providing people with more varied types of limb differences access to this advanced multi-articulating hand.

The COVVI Hand - Flex Wrist

The original 'Flex Wrist' version of the COVVI Hand enables users to adjust their wrist angle for various tasks, offering 27° of flexion and extension that can be securely locked. This allows users to flex inwards to perform tasks like zipping up a jacket, as well as extending outward for activities like pushing carts.



The COVVI Hand - Straight Wrist

The 'Straight Wrist' version of the COVVI Hand represents an industry standard, with its static wrist design.

At approximately 7% lighter than the COVVI Hand 'Flex Wrist' version, installing this variant provides a solution that is both lighter and shorter, while still featuring an Electronic Quick Disconnect (EQD) connection. This model can easily be swapped with other prosthesis tools without clinical assistance.



4.0 The COVVI Hand

The COVVI Hand - Friction Wrist

The ‘Friction Wrist’ version has been designed specifically for people with a longer residual limb.

As the shortest and lightest variant of the COVVI Hand when installed on a prosthesis, its connection can be built directly into a user’s socket. This means that for users with longer limbs, the shorter wrist configuration provides an improved cosmetic appearance and positioning.

In addition to this, the Friction Wrist variant is just under 12% lighter compared to the COVVI Hand when a Flex Wrist is installed.



Structural Height	Small Hand	Medium Hand	Large Hand
COVVI Hand - Flex Wrist	122 + 25mm / 4.8 + 0.98 inches	133 + 25mm / 5.2 + 0.98 inches	136 + 25mm / 5.4 + 0.98 inches
COVVI Hand - Straight Wrist	117 + 25mm / 4.6 + 0.98 inches	128 + 25mm / 5.0 + 0.98 inches	131 + 25mm / 5.2 + 0.98 inches
COVVI Hand - Friction Wrist	117 + 10mm / 4.6 + 0.4 inches	128 + 10mm / 5.0 + 0.4 inches	131 + 10mm / 5.2 + 0.4 inches

The COVVI Hand structural height is defined as the distance from the Thumb Tip to the Base of the EQD Lamination Ring or Friction Wrist Lamination Ring.

The first number in each measurement above is the distance from the Thumb Tip to the Base of the Hand. The second number is the distance from the Base of the Hand to the Base of the Lamination Ring. The overall structural height is the total of these two measurements.


4.0 The COVVI Hand

4.5 Technical Information

COVVI Hand Voltage	7V to 8.4V DC
Maximum Current Consumption	5A (Max. Constant) / 7A (Max. peak, 50ms duration)
Device Weight* COVVI Hand with Flex Wrist	Small: 647g / 1.43 lbs Medium: 672g / 1.48 lbs Large: 703g / 1.55 lbs
Device Weight* COVVI Hand with Straight Wrist	Small: 598g / 1.32 lbs Medium: 624g / 1.38 lbs Large: 655g / 1.44 lbs
Device Weight* COVVI Hand with Friction Wrist	Small: 570g / 1.26 lbs Medium: 593g / 1.31 lbs Large: 620g / 1.37 lbs
Maximum Loading Capacity	90kg / 198 lbs across knuckle
Maximum Carrying Capacity	30kg / 66 lbs across the fingers
Maximum Finger Load Limit	16kg / 35 lbs
Full Open to Full Close to Achieve Grip	Tripod Grip: 0.4 seconds Key Grip: 0.6 seconds Power Grip: 0.8 seconds
Grip Force	Tripod Grip: 35N / 8 lbf Key Grip: 22N / 5 lbf Power Grip: 90N / 20 lbf



Maximum Loading Capacity of push-down force
90 kg / 198 lbs over the knuckles

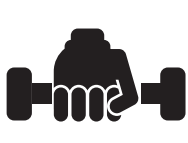

Maximum Carrying Capacity picking weight up
30 kg / 66 lbs



Avoid pressure on thumb

***Device Weight can vary +/- 5 grams**
All conversions into lbs are rounded to two decimal places

5.0 Fitting the COVVI Hand

5.1 Connecting/Disconnecting the COVVI Hand

These instructions only apply to the COVVI Hand if a Flex or Straight Wrist is installed

- Make sure the Power Supply is OFF before connecting/disconnecting the hand to and from the prosthesis.
- Only turn the power on once the hand is properly connected to the prosthesis.

Please see the image below for instructions on how to connect the hand to the socket.

Remove protective cap

Ensure no debris is inside the EQD core before connecting the hand to the socket.




**If the ball bearings are locked out radially, the hand will not plug in.
Please Note: The Wrist Plate marking may vary between variants.**

5.0 Fitting the COVVI Hand



If the COVVI Hand is not used for a long period, protect the EQD from debris with the cap. Do NOT use the COVVI Hand if it is visibly damaged.

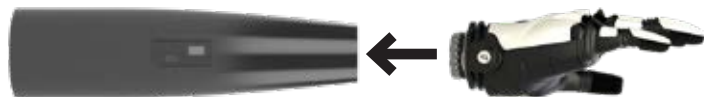

- To connect the hand to the socket, push and twist the EQD into the wrist on the socket.
- To check that the hand is properly connected, twist it in the socket and pull lightly. If it wobbles, it isn't properly connected.



- To disconnect the hand from the socket, twist and pull it.
- If the hand is not properly connected, you may experience BLE connectivity issues or App connectivity issues.

For fitment and removal of the COVVI Hand - Friction Wrist, please refer to the COVVI Hand - Friction Wrist Lamination and Fitment Guide:


<https://www.covvi.com/wp-content/uploads/documents/Friction-Wrist-Lamination-and-Fitment-Guide.pdf>

5.0 Fitting the COVVI Hand

5.2 E-Paper Screen

The COVVI Hand has an E-Paper screen which shows important hand information. It requires no power to display and little power to change the display, so it doesn't drain the battery life.

When you turn the hand on, the E-Paper screen will display the following:

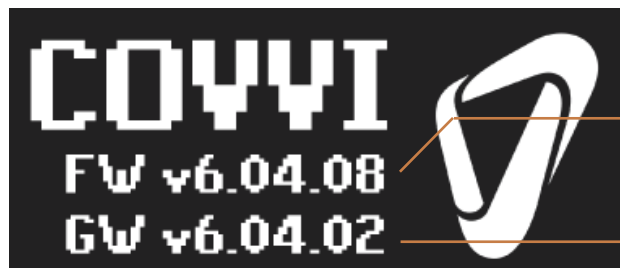

FW: The high-level Firmware. The image shows Version 6 followed by the date it was created.

GW: The low-level Gateway. The image shows Version 6 followed by the date it was created.

It will then display the following:



It always shows the current grip, table and battery life. The Bluetooth icon will be displayed for 60 seconds, and then it will disappear as Bluetooth will automatically turn off; this is a battery saving feature. If you connect to the hand through the App, the Bluetooth icon will remain on the E-Paper screen until you close the App and disconnect from the hand.

5.0 Fitting the COVVI Hand

5.3 System Notifications



When Auto Grip is enabled, if the hand senses an object is slipping it reapplies a small continuous force to grasp the object. Auto Grip works in Power Grip, Tripod Grip, Precision Open and Precision Closed.

Auto Relax moves the hand to a relaxed position after a period of inactivity. A user input signal will cause the hand to revert to its previous grip, and the icon will be removed.

A short press to the dorsal button (under two seconds) puts the hand into Standby Mode, where user inputs from the electrodes have no effect on the hand. Pressing the dorsal button again will exit Standby Mode, and the icon will be removed.

The CAN bus icon is displayed whenever the hand is connected to a CAN-enabled device, such as the COAPT Gen2 system.

Please Note: This icon is only displayed once Bluetooth has turned off (the two icons share the same location on the display). The warning icon shows that the hand has detected a fault. The screen will also display text to identify what the fault is, for example, 1 of 1 Fault Thumb Sensor.

5.0 Fitting the COVVI Hand



When the COVVI Hand is too hot, outside of its normal operating temperature (see section 3.0), the following warning message and icon will appear. If this message appears, please disconnect the hand and allow it to cool down.



When the COVVI Hand is too humid, outside of its normal operating humidity (see section 3.0), the following warning message and icon will appear. If this message appears, please disconnect the hand and allow it to dry in a warm, dry environment.

Any function that refreshes the screen, for example, a grip change, will change the message displayed, but the warning icon will remain visible.

If the warning icons and messages do not disappear, please get in contact with COVVI Customer Services.



The following notifications will appear when you select to update the firmware in COVVI Go. The update downloads from the App to the hand, and then it updates the hand.



5.4 Glove Removal

For the removal and fitting of the standard COVVI Glove, refer to the COVVI Glove Replacement Guide in the User Guides section of the COVVI website: www.covvi.com/user-guides

Please refer to the standard Glove Fitment video for the COVVI Hand with the Flex Wrist, and refer to the Retention Ring Fitment video if your COVVI Hand is fitted with either the Straight or Friction Wrist. Glove fitment between hand variants is the same; however, the glove retention differs between models.

6.0 Grip Patterns

6.1 Opposed Grips

You can find more opposed grips in the User Grips section of the App.



Power Grip



**Precision
Closed Grip**



**Precision
Open Grip**



Trigger Grip



Tripod Grip



Glove Grip



Rock Grip

6.0 Grip Patterns

6.2 Non-Opposed Grips

You can find more non-opposed grips in the User Grips section of the App.



Finger Point Grip



Tap Grip



Relaxed Grip



Mouse Grip



Column Grip



Key/Card Grip




Phone Grip

7.0 The COVVI Power Supply

There are two Power Supply systems available: COVPS-2600 cylindrical cells with a system weight of 107 grams ideal for users with a short residual limb, and COVPS-1600 prismatic flat cells with a system weight of 74 grams and a web connecting the cells. This web allows the batteries to be folded, without bending the cells, to follow the contours of the inside of the prosthesis socket. This system is ideal for users with a long residual limb.

Both Power Supplies come with a Control Panel, which weighs an additional 18 grams. It has USB-C charging, the COVPS-1600 Power Supply takes under three hours to charge and the COVPS-2600 Power Supply takes roughly four hours to charge.

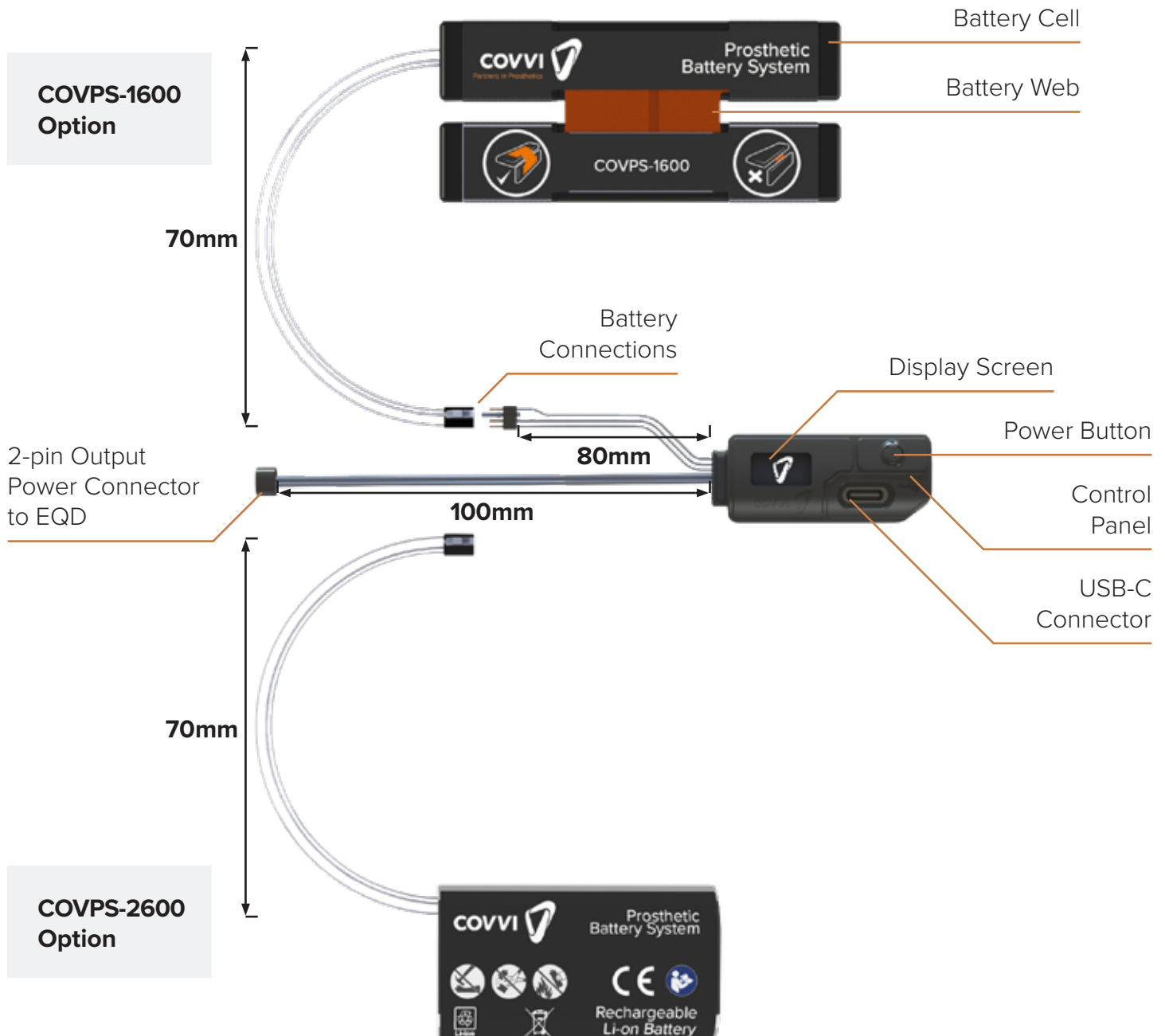


*The Power Supply is optimised for right handed socket installation. For left-handed sockets, the wire must be routed under the Control Panel, leaving 56mm of free wire.

7.0 The COVVI Power Supply

7.1 What's Inside



Cutting Template
(Located in the COVPS Box)

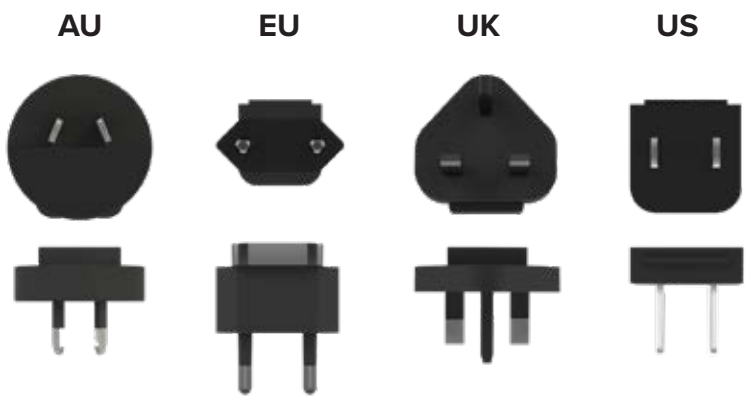

USB-C Wall Charger



Only use with COVVI Batteries

7.0 The COVVI Power Supply

7.2 Safety Precautions

Please read the following safety precautions prior to fitting the COVVI Power Supply.

- The Power Supply should only be fitted by a Certified Prosthetist or Technician.
- This product is not designed to be disassembled or serviced. COVVI Ltd. and COVVI USA INC have the right to void the warranty of all products that have any type of modification or damage caused by any unauthorised or untrained personnel. Any damage caused by intentional harm or neglect will not be covered under the warranty.
- Please make sure the Power Supply is **OFF** before connecting/disconnecting the prosthesis device, as this will avoid any unpredictable spikes in the electric current from the batteries.
- Do **NOT** attempt to use the prosthesis device while the batteries are charging. When the batteries are charging, the power will automatically turn off. If, for any reason, the power does not turn off while charging, using the prosthesis device can be potentially unsafe.
- Do **NOT** use the Power Supply if there is any visible sign of damage to the Power Supply charger, Power Plug and/or Cables.
- Do **NOT** expose the Power Supply to an open flame or submerge it in water. This could damage the screen and affect the battery's ability to hold a charge.
- Do **NOT** use any solvents or abrasives to clean the Control Panel, as this might damage the paint finish, the clarity of the screen or the integrity of the silicone moulding.
- Do **NOT** fit the batteries into a small, enclosed space. The batteries are designed to inflate when they fail. Stopping this inflation may cause the generation of excessive heat, injury, or death.
- Individuals who are exposed to hazardous environments that contain flammable liquid or gas should **NOT** use the COVVI Power Supply when in those environments.
- Ensure access to the Wall Plug to enable easy isolation if required.

Please Note: We recommend you only use the COVVI USB-C Charger to charge the COVVI Power Supply.



This product uses semiconductors that can be damaged by electrostatic discharge (ESD).

7.0 The COVVI Power Supply

7.3 COVPS-1600 Batteries

The COVPS-1600 Battery is made up of four prismatic cells that can be bent either way along the brown connecting web.



7.4 COVPS-2600 Batteries

The COVPS-2600 Battery is made up of two cylindrical cells.



7.0 The COVVI Power Supply

7.5 Activating the System

Press and hold the Power Button on the Control Panel for two seconds to turn the Power Supply on. The screen will display COVVI for two seconds and then the current charge level for five seconds. The system is now providing power to the hand. To turn the Power Supply off, press and hold the Power Button for two seconds.



To reduce power consumption, the screen is blank. To see the charge level, press the Power Button.



7.0 The COVVI Power Supply

7.6 Charge Indication



10 blocks equals 100%



9 blocks greater than or equal to 90%



8 blocks greater than or equal to 80%



7 blocks greater than or equal to 70%



6 blocks greater than or equal to 60%



5 blocks greater than or equal to 50%



4 blocks greater than or equal to 40%



3 blocks greater than or equal to 30%



2 blocks greater than or equal to 20%



1 block greater than or equal to 10%



No blocks less than 10%



7.51V

23°C

When the system is charging, it will cycle through these screens



If there is a fault with the Power Supply, the OLED screen on the Control Panel will show the following warning icon. Please turn off the Power Supply and contact COVVI Customer Services.

7.0 The COVVI Power Supply

7.7 Fitting the COVVI Power Supply

The Control Panel should only be fitted by a Certified Prosthetist or Technician.


The information below is to be used as a guide; individual prostheses will differ from patient to patient.

We recommend that the Control Panel is placed on the upper side of the prosthesis socket, away from any area which may be loadbearing when carrying heavy items. Lamination dummies are provided if recesses are required in the outer socket wall.

Control Panel Aperture



Apply the cutting template to the top side of the prosthesis socket.

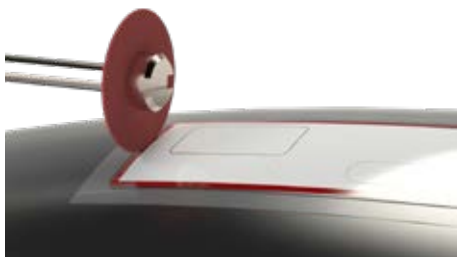

Carefully cut along the red outline with a cutting disc.

File down any sharp edges.



Remove the cutting template and apply a small bead of silicone adhesive to the edge of the aperture.



Insert the wires through the aperture before gently pressing the Control Panel into place.

The batteries can either be connected to the Control Panel prior to being fitted to the limb or after it is fitted if there is sufficient access space.

To ensure correct function, it is important to inspect the Control Panel carefully and ensure that the USB-C port is free of debris. If any debris is discovered, please carefully remove it with a thin, non-metallic object to ensure the port is not damaged.


If the Control Panel needs cleaning, mild soap and a soft, damp cloth should be used. Abrasive cleaners and cleaning cloths will scratch the covers.

7.0 The COVVI Power Supply

7.8 System Specification

Power Supply	COVPS-1600
Typical Battery Capacity	1600mAh
Nominal Voltage	7.4V
Maximum Current Draw	7A Peak
Battery Unfolded / Folded (Inches)	4.1 x 2.0 x 0.4 / 4.1 x 0.8 x 0.8
Battery Unfolded / Folded (mm)	104 x 50 x 9.5 / 104 x 20 x 19.5
System Weight	74 grams (+18 grams for Charge Point)
Humidity	Maximum 80% humidity, non-condensing
Charge Temperature Range	10°C to +45°C
Discharge Temperature Range	-20°C to +60°C

COVPS-1600 (Unfolded)

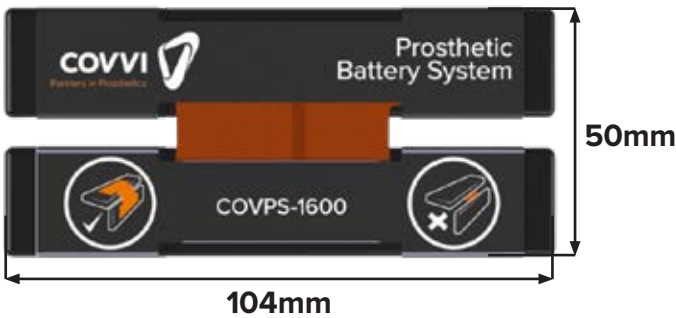

COVPS-1600 (Folded)



Control Panel




Please Note: The battery size in mm is subject to +/- 2mm due to manufacturing tolerances.

7.0 The COVVI Power Supply

Power Supply	COVPS-2600
Typical Battery Capacity	2600mAh
Nominal Voltage	7.4V
Maximum Current Draw	8A Peak
Battery (Inches)	2.8 x 1.6 x 0.8
Battery (mm)	71 x 40 x 20
System Weight	107 grams (+18 grams for Charge Point)
Humidity	Maximum 80% humidity, non-condensing
Charge Temperature Range	10°C to +45°C
Discharge Temperature Range	-20°C to +60°C

COVPS-2600

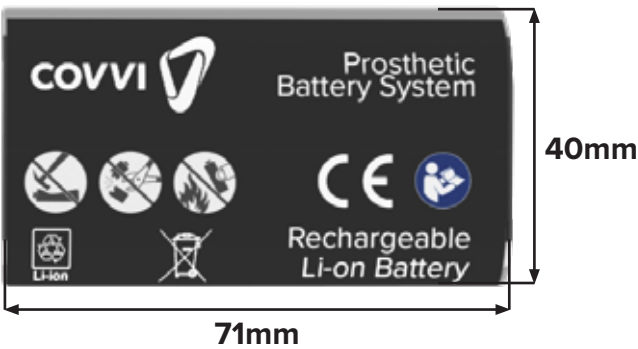

Control Panel



Please Note: The battery size in mm is subject to +/- 2mm due to manufacturing tolerances.

7.0 The COVVI Power Supply

7.9 Battery Replacement

As with all batteries, the capacity will degrade over time; we expect this to become noticeable after 12 months of use. The COVVI Power Supply will support 300 charge and discharge cycles. After this, the time it takes to fully flatten the batteries will become shorter. Replacement cells are available for both Power Supply versions.

It is imperative that you replace the batteries with those of the same capacity. COVSB-1600 is only compatible with the COVPS-1600 Power Supply, and COVSB-2600 is only compatible with the COVPS-2600 Power Supply. Deplete the old batteries until flat, and the limb no longer moves.

Disassemble the limb and replace the batteries. Reassemble the limb, connect the batteries to the USB-C charger and allow them to fully charge. If this initial charge is not conducted properly, it can negatively affect how the system calculates the remaining charge.

Check your local waste regulations to safely dispose of your old batteries.

8.0 Electromagnetic Compatibility (EMC)

Electromagnetic Compatibility, or EMC, refers to the ability of a device to function properly in its electromagnetic (EM) environment without causing or experiencing interference. The COVVI Hand is designed to ensure that it does not emit levels of electromagnetic energy that would cause interference with other nearby electronic devices, nor does it allow electromagnetic interference (EMI) from other devices to affect its performance. International standards and regulatory guidelines help minimise the risk of interference between electronic medical devices.

The COVVI Hand complies with the IEC 60601-1-2 standard, which outlines the requirements for medical electrical equipment in terms of basic safety and essential performance in relation to electromagnetic disturbances. This standard ensures that the device meets the necessary criteria for electromagnetic emissions and immunity to interference.

Please Note: The emission characteristics of the COVVI Hand make it suitable for use in both professional healthcare environments and residential settings, according to CISPR 11 Class B standards. This ensures that the device provides adequate protection against interference with radio communication services. In the unlikely event that interference does occur with nearby radio communication equipment, the user may need to take simple corrective actions, such as relocating or reorienting the device.



To prevent interference or operational issues, avoid placing this device directly next to or stacked with other electrical equipment. If such an arrangement is necessary, ensure both this device and the other equipment are monitored closely to verify that they are functioning correctly and are not affected by interference.



Only use approved accessories, electrodes, and cables provided or recommended by COVVI. Use of non-approved components may increase electromagnetic emissions or reduce the device's electromagnetic immunity, potentially leading to improper operation or device malfunction.



Keep portable RF transmitters (including radios and mobile phones) at least 30cm (12 inches) away from the COVVI Hand. Failure to maintain this distance could result in performance degradation due to interference with the device's electromagnetic sensitivity.

9.0 Component Compatibility



The use of incompatible prosthetic components can cause injury due to unexpected errors. See the table below to show product compatibility. Any device or component that is not listed below is NOT recommended to be used with the COVVI Hand.

Company	Product	Part Number	Compatibility	Recommendations
Electrodes				
Ottobock	Electrode	13E200=50	Yes	N/A
Ottobock	Electrode	13E200=60	Yes	N/A
Ottobock	Suction Socket Electrode	13E202=50	Yes	N/A
Ottobock	Suction Socket Electrode	13E202=60	Yes	N/A
Steeper	Electrode	Elec 50	Yes	N/A
Steeper	Electrode	Elec 60	Yes	N/A
Steeper	Seal-in Electrode	ELSK50	Yes	N/A
Steeper	Seal-in Electrode	ELSK60	Yes	N/A
Ossur	Compact Electrode Kit, 50Hz (300mm Cable)	PL091050	Yes	N/A
Ossur	Compact Electrode Kit, 50Hz (600mm Cable)	PL091127	Yes	N/A
Ossur	Compact Electrode Kit, 60Hz (300mm Cable)	PL091060	Yes	N/A
Ossur	Compact Electrode Kit, 60Hz (600mm Cable)	PL091128	Yes	N/A
IBT	Glide	TBC	Yes	N/A
Batteries				
Ottobock	MyoEnergy Integral	757B35=5	Yes	N/A
Steeper	S-Charge System	SCBP2200	Yes	N/A
Touch Bionics	Replaceable Battery Assembly Kit	PL238149	Yes	N/A
Touch Bionics	Replaceable Battery Assembly Kit w/ switch block	PL238163	Yes	N/A
Ossur	1300mAh battery	PL000336	Yes	N/A
Ossur	2000mAh battery	PL000335	Yes	N/A
IBT	FlexCell	1027201/1027202/ 1027203/1027204	Yes	N/A
Vincent	Vincentaccu Flex	flex1290	Yes	N/A

9.0 Component Compatibility

Company	Product	Part Number	Compatibility	Recommendations
Elbow				
Fillauer	Utah Arm 3 Base	5010110/5010111/ 5010112	Yes	Please contact a COVVI representative for information on compatibility with the Friction Wrist
Fillauer	Utah Arm 3	5010035/5010036/ 5010038	Yes	
Fillauer	Utah Arm 3+	5010039/5010040/ 5010041	Yes	
Fillauer	Utah Hybrid Arm	5010042/5010043/ 5010044	Yes	
Steeper	Espire Pro Elbow	EEP	Yes	
Steeper	Espire Hybrid Elbow	EEH	Yes	
Steeper	Espire Classic Elbow	TBC	Yes	
Steeper	Espire Classic Plus Elbow	EECP	Yes	
Motion Control	Motion Arm EL	5010142/5010144/ 5010143	Yes	
Motion Control	Motion Arm ML	5010139/5010141/ 5010140	Yes	
Ottobock	ErgoArm Hybrid Plus	12K44 (All Sizes)	Yes	
Ottobock	ErgoArm Electronic Plus	12K50 (All Sizes)	Yes	
Ottobock	Analog Adapter	13E100	Yes	Compatible with hands fitted with an EQD, not a Friction Wrist Also needs an analogue adapter
Glaze	Whizzlink		Yes	
Ottobock	DynamicArm Elbow	12K100N (All Sizes)	Yes	
Ottobock	DynamicArm Plus Elbow	12K110N (All Sizes)	Yes	
Connection Cables				
Ottobock	Electrode Cable	13E129 (All Sizes)	Yes	N/A
Ossur	3-way Cable 300mm	PL091029	Yes	N/A
Ossur	3-way Cable 600mm	PL091030	Yes	N/A
Wrist				
Ottobock	MyoRotronic	13E205	Yes	Compatible with hands fitted with an EQD, not a Friction Wrist
Ottobock	Electric Wrist Rotator	10S17	Yes	
Fillauer	MC Standard Wrist Rotator	5010045	Yes	
Fillauer	MC ProWrist Rotator	5010056	Yes	
Fillauer	MC Powered Flexion Wrist	3010993	No	

9.0 Component Compatibility

Company	Product	Part Number	Compatibility	Recommendations
Pattern Recognition				
COAPT	COAPT	G2SC	Yes	N/A
Ottobock	MyoPlus TR	13E520	No	N/A
Lamination Rings				
Ottobock	Lamination Ring	10S1 (All Sizes)	Yes	Compatible with hands fitted with an EQD, not a Friction Wrist
Steeper	Lamination Ring	QDALR-40/45/50	Yes	
Ossur	Lamination Ring	089003/091037	Yes	
Ottobock	Lamination Ring	9S110 (All Sizes)	Yes	Compatible with hands fitted with a Friction Wrist, not an EQD
Coaxial Plug				
Ottobock	Coaxial Plug	9E169	Yes	Compatible with hands fitted with an EQD, not a Friction Wrist
Steeper	Coaxial Plug Kit	QDACO	Yes	
Ossur	4-Way Coaxial Plug	PL091036	Yes	
Coupling Piece				
Ottobock	Coupling Piece	10S4	Yes	Compatible with hands fitted with an EQD, not a Friction Wrist
Steeper	Coupling Piece Kit	QDACP	Yes	
Ossur	Coupling Piece	PL091032	Yes	

All systems which are compatible with a standard EQD. Compatible with a 7-8.4v battery supply. Compatible with 0-5v electrode signals.

WARNING: Use of the COVVI Hand adjacent to or stacked with other equipment not listed should be avoided because it could result in improper operation. If such use is necessary, the COVVI Hand and the other equipment should be observed to verify that they are operating normally.

WARNING: Use of accessories, transducers, and cables other than those specified or provided by COVVI Ltd. and COVVI USA INC could result in increased electromagnetic emissions or decreased electromagnetic immunity of the COVVI Hand and result in improper operation.


WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm (12 inches) to any part of the COVVI Hand, including cables specified by COVVI Ltd. and COVVI USA INC, as this could result in degradation of the performance of this equipment.

10.0 COVVI Go

10.1 COVVI Go

Everything you need to configure a COVVI Hand is within COVVI Go, which is compatible with iOS version 12.0+ and Android version 10.0+. The App has been designed to support you and your patient through every step of your journey with a COVVI Hand. It is an intuitive platform, making programming easy and logical, ensuring you can configure your patient's hand to best suit their daily activities and lifestyle to ensure they get the most out of it. The COVVI Hand comes with a default setup allowing for a smooth adaptation period for new users and efficient interchangeability from other prosthesis devices for existing users.

To download COVVI Go, go to our website www.covvi.com/covvi-go/ to find links to the Apple Store (for iOS) and Google Play Store (for Android), where you can download the App. You can also search for the App directly in both stores by typing in 'COVVI Go'. When the App has downloaded and installed, the icon will appear on your device's home screen.

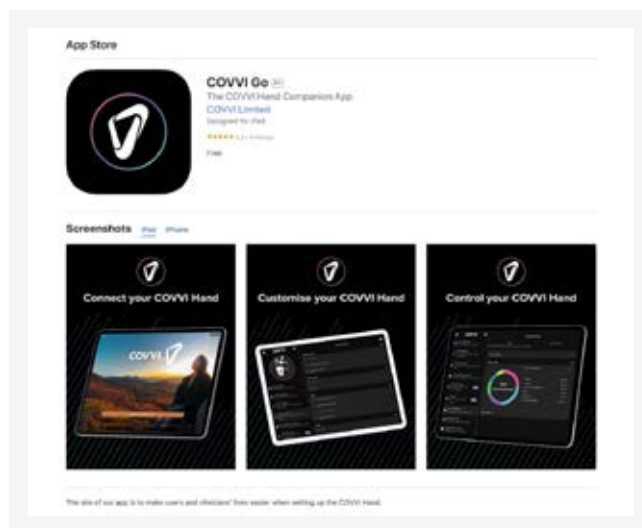

10.2 Registration

You can register within the App by opening it once it has finished installing, clicking on 'Get Started', then on 'Sign Up'. Choose either 'Prosthetic User' or 'Clinician' and then click on 'Next' to fill out all the fields that are relevant to your account type. Once you complete all the required sections, you will receive an email with a link to validate your email and set your password.

10.3 COVVI Portal

You are responsible for adding your patients. To add a new user, go to the COVVI Portal, portal.covvi.com, click on 'New User', fill out their details and click on 'Submit'. There is an option to activate a 'Restricted Mode' to simplify the process for new users, whereby certain features in the App will be disabled. This includes preventing users from accessing some of the more advanced settings within the App, while still allowing them to test out some of the options in the Hand Controls and Input Settings sections. This can also be deactivated at any point.

There is also a 'Virtual Hand' mode, which allows you to trial the App without connecting to a hand.

10.0 COVVI Go

10.4 Hand Selection

When you open the App, click on 'Get Started' and 'Sign In'. There is an option to reset your password in case you forget it.

Once you have logged in, click on 'Connect Hand'. The App will automatically begin searching for your hand. Click on 'Connect' next to the serial number to connect to your hand. The hand communicates using Bluetooth, which only stays activated for 60 seconds as a battery-saving feature. If you fail to connect before 60 seconds, turn the power to the hand off and back on again.

If you experience Bluetooth connectivity issues, hold the Dorsal Button on the back of the hand for 10 seconds. The E-Paper screen will say 'RESET' and then 'HELLO'.

If connectivity issues persist, please ensure the hand doesn't have a low battery and that you don't have Bluetooth interference.



10.5 Input Settings

We advise that the first section of the configuration completed with a user is the Input Settings section, to ensure the user can comfortably open and close the hand before trying out different Grips and Grip Triggers.

A user's signals can be viewed on either the dials or on a graph view. In the 'Input Dials' tab, there is an Input Commands section that will display the triggers that are activated. There are three adjustable thresholds for both Input A and Input B: the 'On', 'Max' and 'Co-Contraction' thresholds.

The 'On' threshold is the signal strength needed for the hand to move. The 'Max' threshold is the signal strength needed for the hand to move at full speed.

The 'Co-Contraction' threshold will only appear if you map Co-Contraction in the Grip Triggers section.

Near the top right of the screen, you can click on the 'Input Settings' tab to see the user's signals on a graph.

Click on the 'Gear' icon in the top right corner to bring up an additional console for more advanced settings. Here you can calibrate the thresholds for both the 'Open' and 'Close' signals. Here, the user has to give the strongest signal they can comfortably produce during the five second timer. Click on 'Confirm' to save the data, or you can reset the test if they're not satisfied with the results.



10.0 COVVI Go

10.6 Grip Tables

There are four tables available, each of which can be populated with up to six grips. All tables can be empty apart from Table A, which must always have one grip set as the hand's primary grip. You can set up the tables with different combinations of grips that suit activities that the user performs in their day-to-day life. For example, you could populate one table with grips used at home, another with grips used at work and another with grips that are helpful for gardening or outdoor work.

Default Grip Tables:

The default setup has two tables populated with the following grips. These are a mix of opposed and non-opposed grips. Tables C & D are empty and therefore disabled.



A	B	C	D
Tripod	Key	Empty	Empty
Power	Finger Point	Empty	Empty
Key	Tripod	Empty	Empty
Finger Point	Power	Empty	Empty
Empty	Empty	Empty	Empty
Empty	Empty	Empty	Empty

Adding Grips:

To add a grip, select the table you would like to populate and click on 'Add New Grip'. Choose the desired grip from the list and then click on the back arrow to return to the table. The grip will now appear in the next available position. Grips are added one by one. If you are unsure of what a particular grip does, select one to see a video showing the product performing the selected grip.

Removing Grips:

To remove a grip, click on it and then select 'No Grip'.

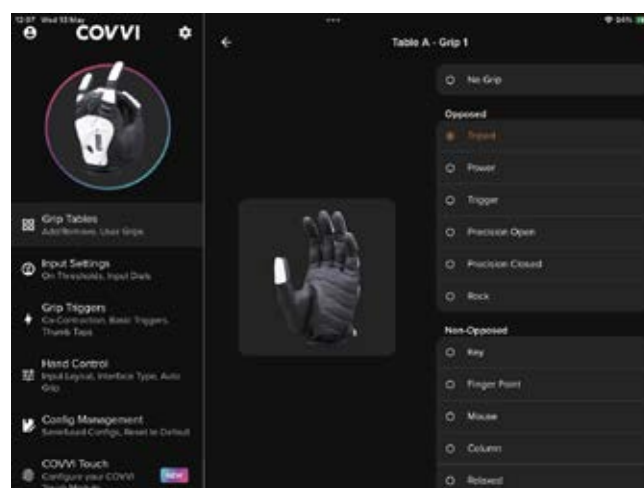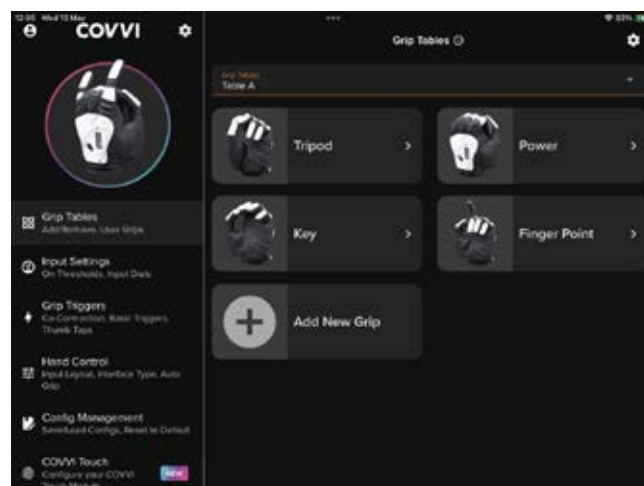


10.0 COVVI Go

10.7 User Grips

There are 14 standard grips available and a library of additional User Grips. To install a User Grip, the hand must have the Firmware Version 4.5.4 or higher. When you add grips to a table, a list will appear divided into different groups. Click on the one you would like to add to your table. You can add six user grips at once.


Click on the back button, and the grip will now appear in the first list of grips under User Grips.

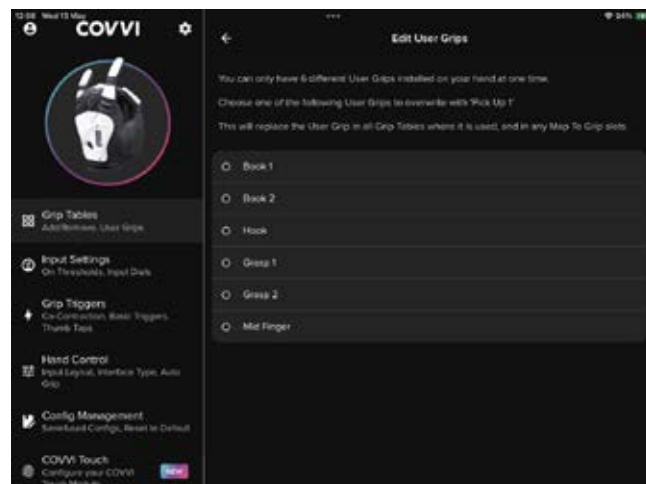

10.8 Grip Triggers

There are 13 available Grip Triggers, 3 myo inputs and 10 passive input triggers on the product, all listed below.

Clicking on any Grip Trigger brings up the list of actions. The action assigned to each Grip Trigger is shown in copper underneath it. You can only assign a single action to a specific Grip Trigger; however, the same action can be set to multiple Grip Triggers if you desire.

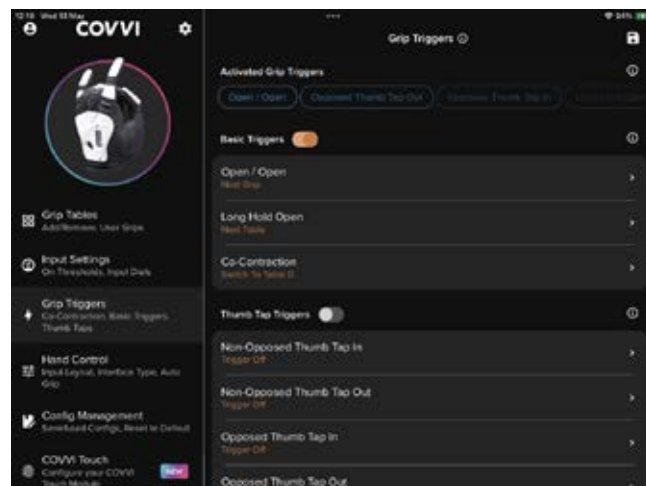

The available Grip Triggers are:

Open/Open - An Open/Open signal is an open signal to fully open the hand, followed by another open signal. The second open signal must be sent before the Open/Open Time.

Long Hold Open - A Long Hold Open signal is an open signal that is maintained. The hand can be closed, partially closed or fully open. The signal must be maintained for the Long Hold Open Time once the hand is fully open.

Co-Contraction - Co-Contraction is when both inputs cross the Co-Contraction thresholds within a certain amount of time from the 'On' threshold being surpassed. The first input to cross the 'On' threshold initiates the Co-Contraction time (a timer customisable by the user), and if both inputs cross their Co-Contraction thresholds before the Co-Contraction time, there will be a successful Co-Contraction trigger.

In the Input Graph, the Co-Contraction thresholds must be higher than the 'On' thresholds. Co-Contraction can be assigned to any action; however, we would recommend that it is mapped to switching between devices, for example, an elbow to a wrist

10.0 COVVI Go

Non-Opposed Thumb Tap In - The thumb is tapped inwards and will move into an opposed position unless assigned to switch to a non-opposed grip, where it will remain in a non-opposed position.

Non-Opposed Thumb Tap Out - The thumb is tapped outwards when in a non-opposed position. If assigned to switch to a non-opposed grip, the thumb will remain in a non-opposed position. If assigned to switch to an opposed grip, the thumb will move into an opposed position.

Opposed Thumb Tap In - The thumb is tapped inwards when in an opposed position. If assigned to switch to an opposed grip, the thumb will remain in an opposed position. If assigned to switch to a non-opposed grip, the thumb will move into a non-opposed position.

Opposed Thumb Tap Out - The thumb is tapped outwards and will move into a non-opposed position unless assigned to switch to an opposed grip, where it will remain in an opposed position.

Thumb Tip Hold - Pressure is applied to the Force Sensitive Resistor in the thumb for two seconds.

Index Finger Tip Hold - Pressure is applied to the Force Sensitive Resistor in the index finger for two seconds.

Middle Finger Tip Hold - Pressure is applied to the Force Sensitive Resistor in the middle finger for two seconds.

Ring Finger Tip Hold - Pressure is applied to the Force Sensitive Resistor in the ring finger for two seconds.

Little Finger Tip Hold - Pressure is applied to the Force Sensitive Resistor in the little finger for two seconds.

Dorsal Button Hold - Hold the dorsal button next to the E-Paper screen for two seconds.

The available actions are:

Trigger Off - No action is assigned.

Next Grip - The next grip in the table is selected.

Previous Grip - The previous grip in the table is selected.

Next Table - The first grip in the next table is selected.

Previous Table - The first grip in the previous table is selected.

Switch to Table A / B / C / D - The first grip in Table A, B, C or D is selected.

Map to Grip - This action allows you to assign the Grip Trigger to go into a specific grip. The grip does not need to be in the grip tables to be mapped.

Toggle Auto Grip - Auto Grip will be turned on and off.

Toggle COVVI Touch Enabled - Enables COVVI Touch vibration feedback function.

Toggle COVVI Touch Repeat - Enables COVVI Touch repeat vibration function.


10.0 COVVI Go

10.9 Hand Control

This section should only need programming once and may not even need any changes at all. The available options are:

- 1. Operating Mode:** One available for dual site, two available for single site, or multiple digital electrodes for use with a pattern recognition system.
- 2. Input Layout:** To order the open and close signals.
- 3. Interface Type:** To choose what control interface is being used.

*If your hand was purchased before April 2024, you may see a different layout.

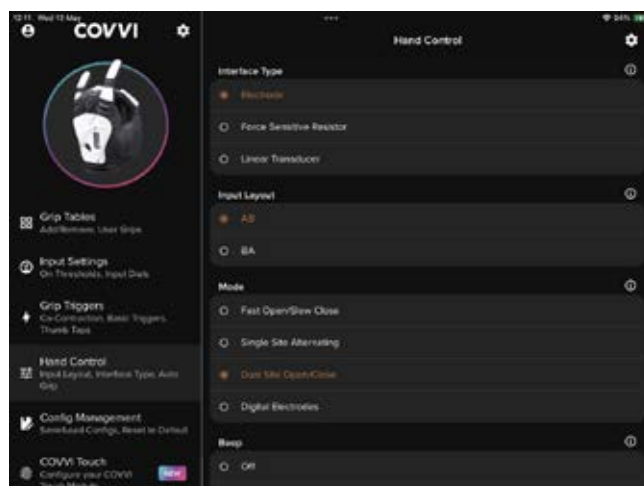

The available options to personalise in the Hand Control section are:

Beep - The beep is an audible alert to notify the user of any activity. It can be assigned to all functions, only power on and standby, only power on or it can be turned off.

Vibration - There is also the option to be alerted to activity with a vibration. Like with the beep, it can be assigned to all functions, only power on and standby, only power on, or it can be turned off.

Grip Settings - Here you can adjust the Grip Strength and Input Speed. The Grip Strength can be proportional to the user's signal or fixed, where it can be manually set to one strength only. The Input Speed is the maximum speed of the digits.

Mode Options - This is dependent on the operating mode being used. You can adjust the time parameters for the Grip Triggers in dual-site use, or the timers for single-site use.

Digit Hold Time - This is the length of time you apply pressure on the FSRs for the Grip Trigger to activate. Moving the slider will increase or decrease the time in 0.5 second intervals, up to a maximum of 5 seconds.

Auto Relax - After a certain period of inactivity, the hand moves into a relaxed position. The time parameter is user-adjustable. Moving the slider will allow users to turn it on and adjust the timer in one-minute intervals, up to a maximum of 10 minutes. Moving it all the way to the left will turn Auto Relax off.

Auto Grip - Auto Grip works in Power Grip, Tripod Grip, Precision Open and Precision Closed. If the fingers sense that an object is starting to slip, they reapply a small continuous force or pulse to gently grasp the object. It monitors the grip every half a second.

Switch to Primary Grip - After a certain period of inactivity, the current grip changes back into the primary grip in Table A. The time parameter is user-adjustable by moving the slider in 15 second intervals, up to a maximum of one minute. Moving it all the way to the left will turn 'Switch to Primary Grip' off.

Thumb Mode - The thumb can be set as one of two options. In 'Latched' mode, the thumb will automatically rotate to the required grip position and stay there for the selected grip. In 'Dynamic' mode, the thumb will rotate to the central position when the hand is open for a more natural look, and then rotate to the required position when forming a grip.

10.0 COVVI Go


Mode Specific Parameters

For Dual Site Open/Close, the following parameters can be established:

Open/Open Time - This is the length of time in which a second open signal must be sent after the hand is fully opened, or after the first open signal, for the Open/Open Grip Trigger to work.

Long Hold Open Time - This is the length of time an open signal must be maintained once the hand is fully opened for the Long Hold Open Grip Trigger to work.

Co-Contraction Time - This is the length of time during which both inputs must pass the Co-Contraction thresholds after the On threshold has been reached.

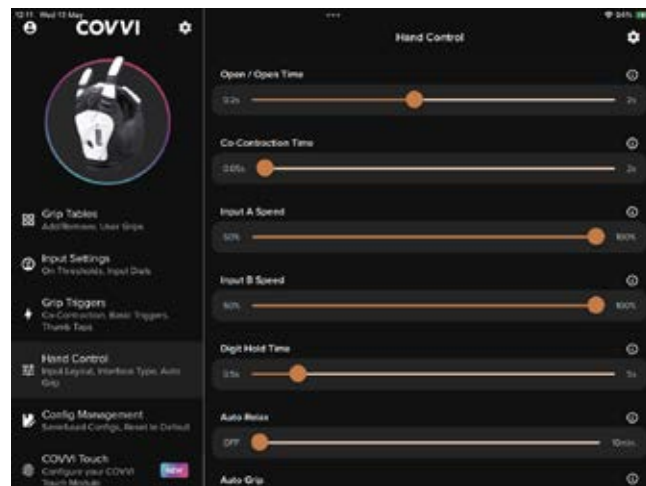

For Fast Open/Slow Close, the following parameters can be established:

Fast Open / Slow Close Time - A signal sent within this time is interpreted as an open signal, and a signal sent after this time is interpreted as a close signal.

Long Hold Open Time - This is the length of time an open signal must be maintained once the hand is fully opened for the Long Hold Open Grip Trigger to work.

For Single Site Alternating, the following parameters can be defined:

Direction Change Time - A signal sent after this time will change the direction of the hand, and a signal sent within this time will keep the direction of the hand the same.

Long Hold Open Time - This is the length of time an open signal must be maintained once the hand is fully opened for the Long Hold Open Grip Trigger to work.

10.0 COVVI Go

10.10 Config Management

In Config Management, you have three options:

1. Save Config allows you to save a new configuration to the cloud or overwrite an existing config, so there is always a backup available. If a config is saved with a name that already exists, a confirmation message will appear to ask if you want it to be overwritten.
2. Default Config Setup allows you to select a COVVI configuration from a range of options, including configurations that are compatible with third-party control systems, such as COAPT.
3. Backup and Restore Live Config is to transfer configurations between devices.



Underneath these options, there is the Config Library, where you can see all configurations that you have saved to the product. To restore a configuration, click on the desired one and then on 'Restore'.

10.11 COVVI Touch

The COVVI Touch module works in tandem with the COVVI Hand to produce vibrational feedback, which is sent to a user's residual limb when their bionic hand comes in contact with an object.

With sensors already present in each digit of the COVVI Hand, COVVI Touch unlocks the technology for a new purpose - letting users know when their bionic hand is touching or holding an object without the need to rely solely on visuals. This provides a better sense of spatial awareness, grasp and even the potential to recognise the size, shape and density of objects.

Settings can be adjusted in this section of the App for elements like deciding which grips to receive feedback in, the strength of the vibration, and how the feedback is relayed (a single pulse, continuous vibration when in contact, or a recurring pulse).




Please Note: COVVI Touch is not compatible when fitted with the COVVI Hand Friction Wrist. It can only be installed on prosthetic systems that feature the COVVI Hand, with either a Flex Wrist or a Straight Wrist.

For more information about being fitted with COVVI Touch, please visit our website: www.covvi.com/covvi-hand/covvi-touch

10.0 COVVI Go

10.12 Games


This section provides a range of training programs to practice different operations of the COVVI Hand. For example, our Training Game is focused on fine-tuning open and close signals by holding a certain signal within a fixed area for different lengths of time. All the games have additional options to customise a user's experience, as well as multiple difficulty options, so users can build up competency with the COVVI Hand at their own pace.

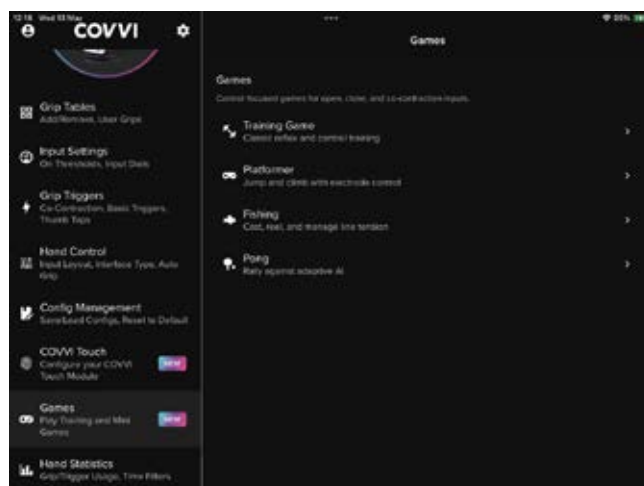

10.13 Hand Statistics

This section displays a range of useful statistics about the COVVI Hand. You can see the top five Grips and Grip Triggers used in the current month, which can also be filtered to show other periods of time as well.


By pressing 'View More', you can see the complete stat history for each individual Grip and Grip Trigger. Selecting any one of these then shows you the detailed history of this information, displayed on a graph.

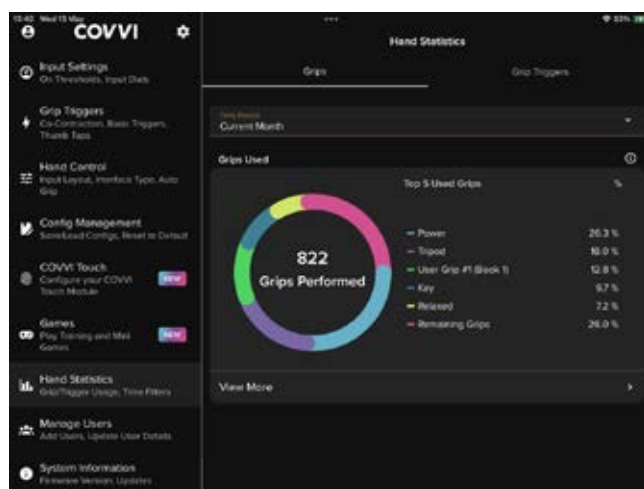

10.14 System Information

In this section of the App, you can find basic information about the COVVI Hand, including the hand's Serial Number and the current Firmware Version.

There is also the option to 'Check For Updates', which will look to see if there are any software updates that are currently available. If there are any software updates, the App will then guide you through this process of updating.

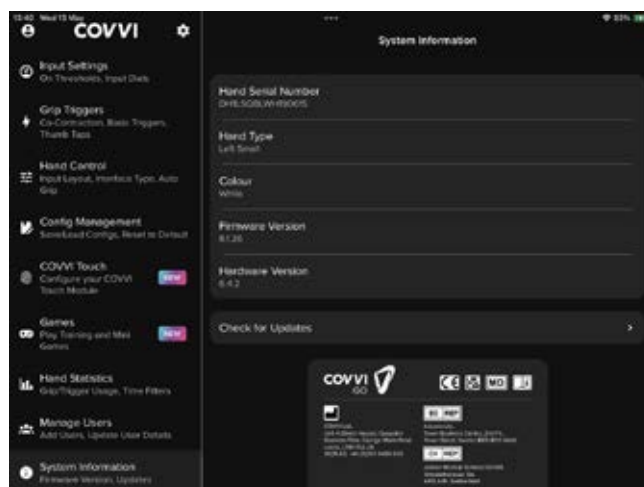

10.0 COVVI Go


10.15 Software Update

To continuously improve the COVVI Hand, an integrated update feature runs on the device firmware. These updates include improvements, new features, bug fixes, etc... These updates do not require the user to send their hand into a COVVI repair centre to be updated; they can be downloaded to the hand remotely from the cloud.

Click on 'Check For Updates'. If there are updates available, you must update the Firmware first, then the Gateware.

To perform the update, click on it and select 'Download and Install'. The update first downloads from the App onto the hand, then the update is installed to the hand. This process can take up to 30

minutes, so we recommend keeping the hand close to the device to ensure a strong Bluetooth connection is maintained, and download speeds are optimal. The E-Paper screen will change to indicate an update is taking place. The hand must remain connected to a power source until the update is complete. If it is disconnected from power, it will stop the update, which means it will need to start again. It cannot be connected to a charger during the update because this will turn the power off to the hand. Once the update is complete, the hand will reboot. Repeat the process for the Gateware.

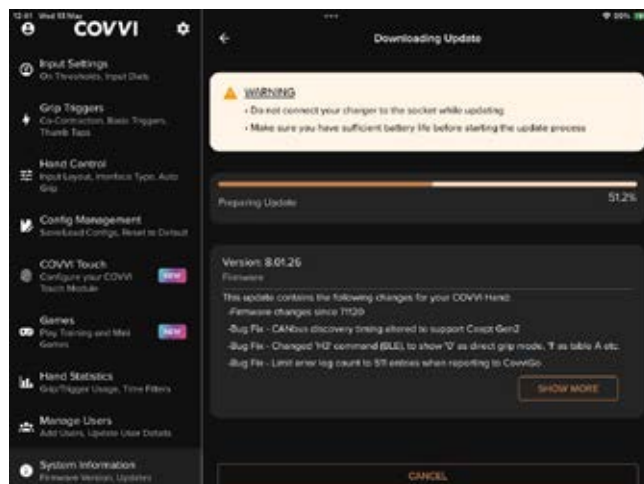


10.16 Remote Assist

The COVVI Remote Assist feature allows a COVVI Representative to provide real-time support to clinicians and users to implement instant configuration changes to the live configuration, software updates, and perform troubleshooting to any COVVI Hand, anywhere in the world. Accredited clinicians have access to this platform to do the same for users.

From the main menu of the App, you can access the Remote Assist feature by tapping the 'Gear' icon in the top right corner, and then selecting 'COVVI Remote Assist' from the Console.

Please Note: We can **ONLY** connect to the hand if a clinician and/or COVVI is notified first, and a meeting time is confirmed. Ensure that the window is kept open in the App while we connect to the hand. You will know when we are connected, as under technician status, it will say 'Connected' in green.

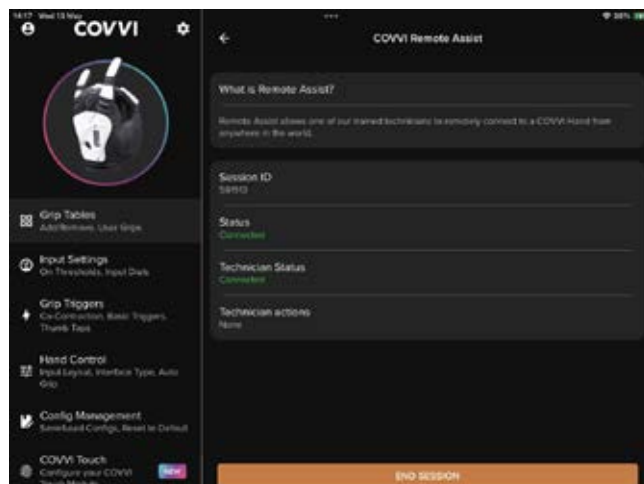

11.0 Warranty

Please read these Warranty Terms and Conditions carefully before using all COVVI products and services. We may modify these Terms for any reason at any time, by posting a new version on our website; these changes do not affect rights and obligations that were defined prior to such changes. There are various limited warranties that apply to manufactured products and goods purchased from COVVI.

11.1 24-Month Manufacturer's Warranty

COVVI offers a 24-month manufacturer's warranty for the hand and a 3-month manufacturer's warranty for the glove, both included with the purchase of the hand, which take effect from the date of invoice. If no Extended Warranty is purchased, all of COVVI's obligations shall cease upon the expiration of the manufacturer's warranty.

11.2 Optional Extended Warranty

COVVI offers an extended warranty for an additional year, covering the hand for 36 months, or an additional three years covering the hand for 60 months. The extended warranties have an additional cost and must be purchased with the hand.

11.3 Additional Information

Products are subject to evaluation for warranty. COVVI is not responsible for normal wear, and/or damage caused by excessive force, and/or excessive usage beyond the technical design and/or beyond its reasonable means. COVVI warrants its products against defects in material and workmanship within the warranty period.

Limitation in those instances where changes, alterations or modifications are made in materials at the request or instruction of the customer, the customer agrees not to claim or commence suit against COVVI based on any such disclaimed warranties.

Our obligation is limited only to the repair and replacement of defective parts within the warranty period or, at the sole discretion of COVVI, to refund the purchase price of a full refund, partial refund, or no refund, depending on the condition of the return. The possible refund will be given, subject to our Quality Department.

12.0 Warranty Service and Repairs

Only the COVVI Hand is serviceable; all COVVI supporting products are designed to be maintenance-free. When any COVVI product is sent in for repairs or to be evaluated for faults, it must be returned to COVVI completely assembled and not tampered with. COVVI has the right to void the warranty of all products that have any type of modifications or damage caused by any unauthorised or untrained personnel. The manufacturer's and extended warranty do not cover repairs needed as a result of:

- Any intentional damage or neglect
- Any alterations or attempted repairs without written consent from the Supplier
- Any damage or defects arising due to not following the safety precautions as stated in the Instructions for Use and Technical Manual
- Any damage or defects arising from improper installation
- Any damage or defects arising from using the Hands with incompatible products
- Any damage or defects arising due to using, transporting or storing the Hands outside of the parameters as stated in the Environmental Conditions section of the Instructions for Use and Technical Manual
- Any damage or defects arising due to any illegal activities involving the Hands
- Any damage or defects arising due to accidents while operating a motor vehicle using the Hands
- Any damage or defects due to use by anyone other than the owner

COVVI reserves the right not to repair any Hands under the Extended Warranty if the Customer makes any further use of such Hands after giving notice of a defect.

Should the product be under warranty and have any manufacturing and/or workmanship defects, the part(s) will be replaced at no charge. The warranty holder will be informed of any and all defects that are covered and not covered under warranty. Any product out of warranty that is sent in to COVVI will be evaluated by our Quality Department and quoted with repair parts and labour as required. Once the repair estimate is sent to the customer, they will have a 60-day approval period. Once the 60-day approval period expires, COVVI has the right to ship the device back to the customer.

12.1 Limited Liability

To the fullest extent permitted by law COVVI Ltd. and its affiliates, directors, officers, employees, partners, contractors, or agents will not be liable for any losses or damages whether direct, indirect, incidental, special, punitive, or consequential resulting from the use of the COVVI Hand and associated products, irrespective of whether the Clinician or User has been advised or otherwise might have anticipated the possibility of such loss or damage.

COVVI Ltd. and its affiliates, directors, officers, employees, partners, contractors, or agents shall not be responsible for strikes, labour slowdowns, war, terrorism, riots, severe weather conditions, natural disasters, acts of God or any other forces beyond the reasonable control of COVVI which may result in direct, indirect, incidental, special, punitive, or consequential losses or damage.

13.0 Device Problem Reporting

In the event of any serious incident involving the product, it is essential that it is reported promptly to the manufacturer and the competent authority of the member state in which the user is established, to maintain the highest standards of safety and compliance. Please report any incident to COVVI Ltd. or COVVI USA INC by emailing customerservice@covvi.com.

Please ensure you include the serial number of the COVVI Hand(s) affected when you reach out to us, which can be found on the base of the device, starting with 'CV' and continuing clockwise. If you have a loaner or demo hand, then this serial number begins with 'LN'.

Ensuring that any issues are reported helps maintain the highest standards of safety and compliance.

Flex Wrist



Straight Wrist



Friction Wrist



13.1 Returns

The standard return period for any COVVI product is 14 days from the date of receipt at the shipping address.

If you or your patient would like to request a refund for any COVVI products, you must complete the Return Form on the COVVI website (www.covvi.com/contact/return-request) explaining the reason, and the items must be shipped in their original condition and packaging to COVVI for inspection.

COVVI is under no obligation to refund:

- Any Goods that have been used by your patient.
- Any Goods that are damaged after delivery to the shipping address.
- Any requests received after the standard return period.

The returns policy only applies to products manufactured by COVVI Ltd. and COVVI USA INC and does not apply to products distributed by COVVI Ltd. and COVVI USA INC.

14.0 Disposal

COVVI products should not be disposed of with regular household waste. Proper disposal of COVVI products will avoid having a detrimental impact on health and the environment. For more details on collection and proper disposal, please check with your local authority.


COVVI packaging should be disposed of through your local facilities. By disposing of the packaging in the proper manner, you help to avoid possible hazards for the environment and public health. Observe national registration requirements for importers according to the EU WEEE Directive 2012/19/EU.

All the components of this product and articles contained within have been chosen and evaluated in accordance with compliance with the EU RoHS directive 2011/65/EU: The products are produced with restricted use of hazardous substances to protect the environment.

All the components of this product and articles contained within have been chosen and evaluated in accordance with compliance for the Regulation (EC) 1907/2006 of the European Parliament, "Registration, Evaluation, and Authorization of Chemicals (REACH): The product and articles contained within **DO NOT CONTAIN** any of the 247 REACH - Substances of Very High Concern (SVHC) as updated by ECHA on January 21, 2025.

For the full declaration of conformity statements, please use the links below to view the documents:

www.covvi.com/about/declaration-of-conformity

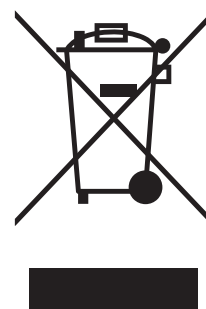

15.0 Symbols Used

Below are the symbols used to convey information about the COVVI product range, which is derived from standards EN 980:2008, ISO 15223, IEC 60417 and IEC 60878:



Indicates that the product conforms with the essential requirements and provisions of MDR 2017/745.



Indicates that the user should refer to the instruction manual or consult the instructions for use before operating the medical device. It directs the user to read the accompanying documentation for safety, operational, or maintenance instructions.



Identifies the entity responsible for the design, production, packaging, and labelling of the device.



A unique identifier assigned to each device that differentiates it from others. Used for tracking and recalls.



A unique identifier assigned to each device that differentiates it from others. Used for tracking and recalls.



Indicates that the device is part of a specific lot or batch, and the accompanying number identifies that particular lot or batch for tracking purposes.



Indicates that the device is a specific product from the manufacturer's catalogue, and the accompanying number identifies that particular product for reordering or reference purposes.



Indicates that the product is classified as a medical device under relevant regulatory requirements, such as EU Medical Device Regulation (MDR) 2017/745. It helps distinguish medical devices from other products.



Indicates that the device must be protected from water or moisture exposure, which could damage the internal components.


15.0 Symbols Used




Indicates that the device should not be disposed of with regular household waste and must comply with WEEE regulations for proper electronic waste disposal.



Specifies the temperature range within which the device can be safely stored or operated.

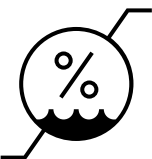

Indicates the range of humidity the device can be safely exposed to or operated within.

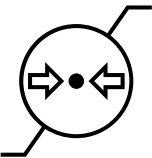

Specifies the range of atmospheric pressure to which the device can be exposed.



Indicates the name and address of the Authorised Representative in the European Union (EU) for the medical device. The authorised representative acts on behalf of the manufacturer in ensuring that the manufacturer meets regulatory requirements in Europe, especially for non-EU manufacturers.



Indicates the name and address of the Authorised Representative in the Switzerland for the medical device. The authorised representative acts on behalf of the manufacturer in ensuring that the manufacturer meets regulatory requirements in Switzerland.



The UKCA marking is the conformity marking used for products being placed on the market in Great Britain.

15.0 Symbols Used




Indicates that the electronic device complies with the FCC regulations regarding electromagnetic interference (EMI) and radio frequency emissions, as defined under Title 47 CFR Part 15 (and other applicable parts)



Indicates that the device has a Type BF (Body Floating) applied part that provides a degree of protection against electric shock.



Non-ionising radiation

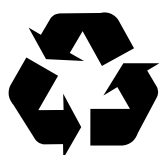

Li-ion

Indicates that the product contains a Lithium-Ion (Li-ion) battery and must be recycled properly. It is a general environmental symbol that emphasise the need to dispose of the battery in an environmentally responsible manner.



This symbol warns against crushing the device or battery. Crushing the device could cause internal damage, including potential battery rupture, which may lead to fire or explosion.



Indicates that the device or battery should not be disassembled. Opening or tampering with the internal components can result in dangerous conditions, such as exposure to chemicals, electric shock, or fire hazards.



Warns against exposing the battery or device to fire or high temperatures. Doing so could cause the battery to overheat, catch fire, or even explode.




A line with a solid line above a dashed line indicates that the device operates on direct current (DC).



A battery icon with a lightning bolt or an arrow indicates that the device or battery is charging or capable of being recharged.

15.0 Symbols Used

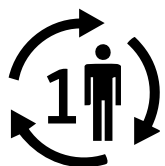

Indicates that the device is intended for use by a single patient, but for multiple times.



Indicates to not load or transport package if damaged.



Indicates that a medical device that should not be used if the package has been damaged or opened and that the user should consult the instructions for use for additional information.



Indicates that the packaging material or product can be collected, processed, and reused in recycling systems rather than disposed of as waste. The symbol encourages proper disposal through recycling channels according to local regulations.



Indicates that the user must consult the instructions for use (IFU) or user manual for important information on safe and effective use, warnings, precautions, or installation and maintenance procedures.



Indicates that the product conforms to the applicable Technical Regulations of Ukraine, as established by the Cabinet of Ministers of Ukraine, and that the manufacturer has successfully completed the required conformity assessment procedure. It confirms that the product meets essential safety and performance requirements and can be lawfully placed on the market of Ukraine.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

- United States/FCC ID: A8TBM71S2
- Canada IC: 12246A-BM71S2
- Europe/CE
- China/SRRC: CMIIT ID: 2016DJ5890



Japan: 005-101150



South Korea: MSIP-CRM-mcp-BM71BLES1FC2



CCAN16LP0011T7

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and **The Electromagnetic Compatibility Regulations 2016 (SI 2016 No. 1091)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Hand
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN001WU
Global Medical Device Nomenclature Code:	41497
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand is intended for individuals with unilateral or bilateral upper limb loss or congenital limb deficiency. It is designed exclusively for use in exoprosthetic fittings of the upper limbs, providing functional and assistive support to users.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Hand Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Hand.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices and Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast)** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Hand
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN001WU
Global Medical Device Nomenclature Code:	41497
European Medical Device Nomenclature Code:	Y061899
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand is intended for individuals with unilateral or bilateral upper limb loss or congenital limb deficiency. It is designed exclusively for use in exoprosthetic fittings of the upper limbs, providing functional and assistive support to users.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Hand Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Hand.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Wrist Lamination Kit (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and **The Electromagnetic Compatibility Regulations 2016 (SI 2016 No. 1091)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	Wrist Lamination Kit
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN003WY
Global Medical Device Nomenclature Code:	41086
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Wrist Lamination Kit is intended to transmit power and control signals to upper limb exoprosthetic terminal devices.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Wrist Lamination Kit Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Wrist-Lamination-Kit.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Wrist Lamination Kit (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices and Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast)** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	Wrist Lamination Kit
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN003WY
Global Medical Device Nomenclature Code:	41086
European Medical Device Nomenclature Code:	Y061899
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Wrist Lamination Kit is intended to transmit power and control signals to upper limb exoprosthetic terminal devices.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-0000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Wrist Lamination Kit Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Wrist-Lamination-Kit.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Power Supply (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and **The Electromagnetic Compatibility Regulations 2016 (SI 2016 No. 1091)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Power Supply
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN004X2
Global Medical Device Nomenclature Code:	36534
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Power Supply is intended to power externally controlled upper limb prosthetic hands.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Power Supply Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Power-Supply.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Power Supply (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices and Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast)** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Power Supply
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN004X2
Global Medical Device Nomenclature Code:	36534
European Medical Device Nomenclature Code:	Y061899
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Power Supply is intended to power externally controlled upper limb prosthetic hands.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-0000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Power Supply Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Power-Supply.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The Glove Kit (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	Glove Kit
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN1GK2Q
Global Medical Device Nomenclature Code:	64736
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Glove Kit is an accessory intended to be worn with the COVVI Hand to protect it from external contaminants and maintain device integrity.
Medical Device Classification:	Class 1 [Rule 1]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full Glove Kit Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Glove-Kit.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The Glove Kit (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	Glove Kit
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN1GK2Q
Global Medical Device Nomenclature Code:	64736
European Medical Device Nomenclature Code:	Y090312
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Glove Kit is an accessory intended to be worn with the COVVI Hand to protect it from external contaminants and maintain device integrity.
Medical Device Classification:	Class 1 [Rule 1]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full Glove Kit Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Glove-Kit.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Go (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Go
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN012WZ
Global Medical Device Nomenclature Code:	64734
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Go is software intended to assist users of the COVVI Hand by enabling the selection and customization of grip functions and the adjustment of control settings.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Go Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Go-App.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Go (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Go
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN012WZ
Global Medical Device Nomenclature Code:	64734
European Medical Device Nomenclature Code:	Z129082
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Go is software intended to assist users of the COVVI Hand by enabling the selection and customization of grip functions and the adjustment of control settings.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Go Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Go-App.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand - Straight Wrist (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and **The Electromagnetic Compatibility Regulations 2016 (SI 2016 No. 1091)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Hand - Straight Wrist
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN027XE
Global Medical Device Nomenclature Code:	41497
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand - Straight Wrist is an external, motor-driven prosthetic hand with a Straight Wrist designed to assist individuals with unilateral or bilateral upper limb deficiencies. Specifically crafted with anatomically correct dimensions, the COVVI Hand - Straight Wrist facilitates activities of daily living through controlled movement of programmed grip patterns.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Hand - Straight Wrist Declaration of Conformity (UK):
www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Hand-Straight-Wrist.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand - Straight Wrist (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices and Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast)** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Hand - Straight Wrist
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN027XE
Global Medical Device Nomenclature Code:	41497
European Medical Device Nomenclature Code:	Y061899
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand - Straight Wrist is an external, motor-driven prosthetic hand with a Straight Wrist designed to assist individuals with unilateral or bilateral upper limb deficiencies. Specifically crafted with anatomically correct dimensions, the COVVI Hand - Straight Wrist facilitates activities of daily living through controlled movement of programmed grip patterns.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-0000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Hand - Straight Wrist Declaration of Conformity (EU):
www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Hand-Straight-Wrist.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand - Friction Wrist (UK, UKCA Marking)

We hereby declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and **The Electromagnetic Compatibility Regulations 2016 (SI 2016 No. 1091)**, and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Hand - Friction Wrist
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN028XG
Global Medical Device Nomenclature Code:	41497
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand - Friction Wrist is an external, motor-driven prosthetic hand with a Friction Wrist designed to assist individuals with unilateral or bilateral upper limb deficiencies. Specifically crafted with anatomically correct dimensions, the COVVI Hand - Friction Wrist facilitates activities of daily living through controlled movement of programmed grip patterns.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Go Declaration of Conformity (UK):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-UK-COVVI-Hand-Friction-Wrist.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

The COVVI Hand - Friction Wrist (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices and Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast)** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Hand - Friction Wrist
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN028XG
Global Medical Device Nomenclature Code:	41497
European Medical Device Nomenclature Code:	Y061899
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Hand - Friction Wrist is an external, motor-driven prosthetic hand with a Friction Wrist designed to assist individuals with unilateral or bilateral upper limb deficiencies. Specifically crafted with anatomically correct dimensions, the COVVI Hand - Friction Wrist facilitates activities of daily living through controlled movement of programmed grip patterns.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-0000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Hand - Friction Wrist Declaration of Conformity (EU):
www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Hand-Friction-Wrist.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:
COVVI Glove Kit 2 (UK, UKCA Marking)

We hereby declare, under the manufacturer’s sole responsibility, that the products described in this document are in conformity with the applicable provisions of the **Medical Devices Regulations 2002 (SI 2002 No. 618, as amended)** and are manufactured and released in accordance with the Quality Management System in force on the date of issue of this declaration. The UKCA marking may be affixed.

Product Name:	COVVI Glove Kit 2
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
EU Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN2GK2V
Global Medical Device Nomenclature Code:	64736
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Glove Kit 2 is an accessory intended to be worn with the COVVI Hand (Friction and Straight Wrist) to protect it from external contaminants and maintain device integrity.
Medical Device Classification:	Class 1 [Rule 12]
Notified Body:	N/A
UKCA Certificate:	N/A
UK Approved Body Name:	Not Applicable for Class I
UK Approved Body Identification Number:	Not Applicable for Class I
Medical Device Regulation Assessment Route:	In conformity with the applicable requirements of the Medical Devices Regulations 2002 (SI 2002 No. 618, as amended), commonly referred to as UK MDR 2002.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Glove Kit 2 Declaration of Conformity (UK):
www.covvi.com/wp-content/uploads/documents/DoC/ Declaration-Of-Conformity-UK-COVVI-Glove-Kit-2.pdf

16.0 Declaration of Conformity

Declaration of Conformity for:

COVVI Glove Kit 2 (EU, CE Marking)

We herewith declare, under the manufacturer's sole responsibility, that the products described in this document are in conformity with the relevant requirements of **Regulation (EU) 2017/745 of the European Parliament and of the Council of 5 April 2017 concerning Medical Devices** and are manufactured and released according to the provision of the quality management system in force on the date of issue of this declaration and the CE Mark may be affixed.

Product Name:	COVVI Glove Kit 2
Legal Manufacturer:	COVVI Ltd. , Unit 4 (Direct House), Quayside Business Park, George Mann Road, Leeds, LS10 1DJ, United Kingdom
Single Registration Number (SRN):	GB-MF-000033926
Basic UDI-DI:	506072613GMN2GK2V
Global Medical Device Nomenclature Code:	64736
European Medical Device Nomenclature Code:	Y090312
Applicable Legislations And Standards:	As per Appendix I - Applicable Legislations and Standards
Configurations And Variants:	As per Appendix II - Product Listing/Schedule
Intended Purpose:	The COVVI Glove Kit 2 is an accessory intended to be worn with the COVVI Hand (Friction and Straight Wrist) to protect it from external contaminants and maintain device integrity.
Medical Device Classification:	Class 1 [Rule 13]
Notified Body:	N/A
CE Certificate:	N/A
EU Authorised Representative:	Advena Limited , Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013, Malta
EU Authorised Representative SRN:	MT-AR-000000234
Medical Device Regulation Assessment Route:	In conformity with Annexes II and III, the Declaration of Conformity has been drawn up in accordance with Article 19 of Regulation (EU) 2017/745.

Full Declaration of Conformity Document

Please click the following link to view the full COVVI Glove Kit 2 Declaration of Conformity (EU):

www.covvi.com/wp-content/uploads/documents/DoC/Declaration-Of-Conformity-EU-COVVI-Glove-Kit-2.pdf



v1.5.2 | June 2026

Uncontrolled Copy when printed
CV-000918-TC, Revision 7



Phone:

Email:

Website:

COVVI Ltd., Unit 4 (Direct House),
Quayside Business Park, Leeds, LS10 1DJ

+44 (0)20 3949 9500

customerservice@covvi.com

www.covvi.com

EU REP

CH REP

Advena Ltd. Tower Business Centre, 2nd Flr.,
Tower Street, Swatar, BKR 4013 Malta

Johner Medical Schweiz GmbH,
Tafelstattstrasse 13a, 6415 Arth, Switzerland